

NovaLife Soft Convex

Features & Benefits

NovaLife 1 Open Soft Convex Clear or with EasiView™ Viewing Option

NOVALIFE	MIDI	MIDI	MIDI	MIDI	MAXI	MAXI	MAXI	MAXI	MAXI	MAXI
	EasiView™	Clear								
Starter hole	REF	REF								
15-24 mm	981-24	982-24			983-24	984-24				
25 mm	981-25	982-25			983-25	984-25				
15-34 mm			981-34	982-34	983-34	984-34				
30 mm			981-30	982-30	983-30	984-30				
35 mm			981-35	982-35	983-35	984-35				
15-44 mm							981-44	982-44		
40 mm							981-40	982-40		
45 mm							981-45	982-45		
15-54 mm									981-54	982-54

Specifications

Convexity depth	6 mm				
Skin barrier	107 x 91 mm	130 x 110 mm	130 x 110 mm	130 x 110 mm	150 x 127 mm
Length x width	239 x 144 mm	239 x 144 mm	259 x 158 mm	259 x 158 mm	259 x 158 mm
Volume capacity	350 ml	360 ml	540 ml	520 ml	500 ml
Box of	10 pcs				

NovaLife 1 Closed Soft Convex With EasiView™ Viewing Option

NOVALIFE	MIDI	MIDI	MAXI	MAXI
	EasiView™	EasiView™	EasiView™	EasiView™
Starter hole	REF	REF	REF	REF
15-24 mm	971-24			
25 mm	971-25			
15-34 mm		971-34		
30 mm		971-30		
35 mm		971-35		
15-44 mm			971-44	
40 mm			971-40	
15-54 mm				971-54

Specifications

Convexity depth	6 mm	6 mm	6 mm	6 mm
Skin barrier	107 x 91 mm	130 x 110 mm	130 x 110 mm	150 x 127 mm
Length x width	239 x 144 mm	239 x 144 mm	259 x 158 mm	259 x 158 mm
Volume capacity	450 ml	430 ml	590 ml	570 ml
Box of	10 pcs	10 pcs	10 pcs	10 pcs

Dansac A/S
Lille Kongevej • DK-3480 Fredensborg • Danmark
Tel +45 4846 5000 • Fax +45 4846 5010
www.dansac.com

NOVALIFE
The next generation

dansac

P04-560P-300 01/2016 © 2016 Dansac A/S onthefront.dk

NEW
Added pouches

NovaLife **Soft Convex**

– For confidence and security

Low headspace

The pouch area above the barrier has been minimised, which reduces the height of the pouch thus making it more discreet and preventing it from slouching.

Release tab

The barrier can be simply and securely removed thanks to the easy-to-see orange non-adhesive release tab.

Oval shaped barrier

The shape of the barrier suits the natural curves of the body and ensures the adhesive area is positioned correctly.

Gentle pressure in the right place

A soft and flexible convex insert enables the appliance to put gentle pressure close to the stoma to help push the stoma into the pouch and to create an enhanced skin seal.




Skin-friendly barrier

The composition of the GX Hydrocolloid barrier is based on natural ingredients and has proven its skin-friendly properties over many years.


Tapered barrier

The surface of the barrier is thickest around the stoma then thins towards the edges to provide secure and gentle adhesion to the skin.

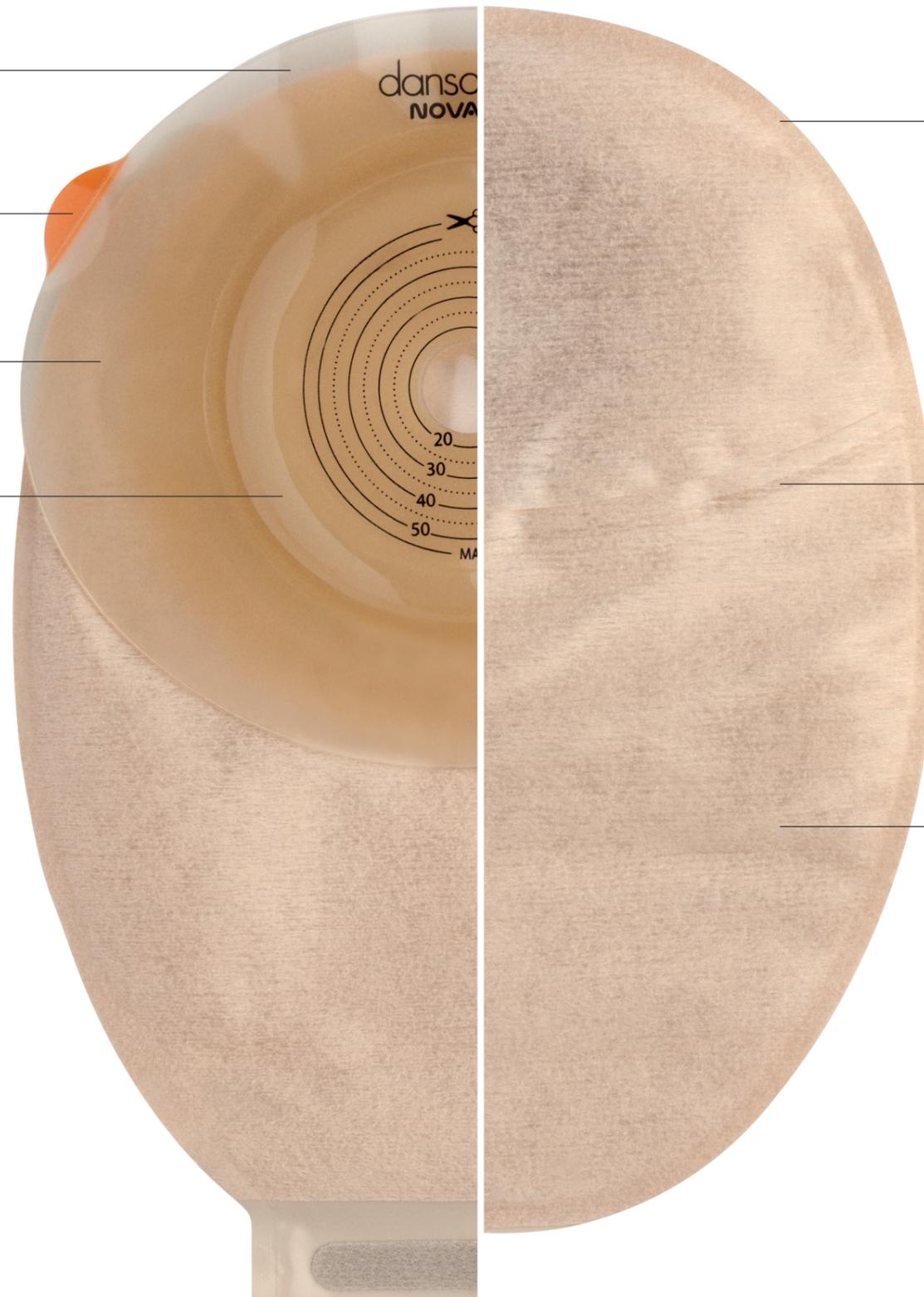

Highly effective NovaLife filter

Membranes on both sides of the NovaLife filter protect it from the contents of the pouch as well as from water outside and make filter patches unnecessary. The filter can expel a large amount of gas very quickly, which minimizes the risk of the pouch ballooning.

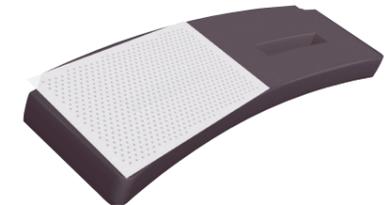

EasiView™ - viewing option



Allows easy inspection of the stoma and facilitates positioning of the pouch.



Water repellent cover

The soft fabric of the pouch is water repellent and will stay dry even when bathing.

Integrated belt ears

Strong, yet flexible belt ears make it easy to attach a Dansac belt to the pouch if more pressure is needed.



NEW
Added pouches

NOVALIFE
The next generation

dansac